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Preface

The staff at JinHao Company would like to thank you for your selection of this product and wishes you good health. With the development of market economics, labour intensity of people is increasing and living rhythm is quickening, so people often feel fatigue and sometimes catch disease. In order that you keep good health , our company will provide you with a electric stimulator.

Our company also offers you a fine service as well.

Thank you!

(Read this manual in detail before using it)


Contents

I . Introduction	17
II . Operating Panel	19
III . Operating Step	20
✧ Operating	20
✧ Attached explanation for the use of remote controller	20
IV . Prescription Functions	20
V . Contraindication	20
VI . Caution	20
VII . Main Technical Parameters	20

1. Introduction

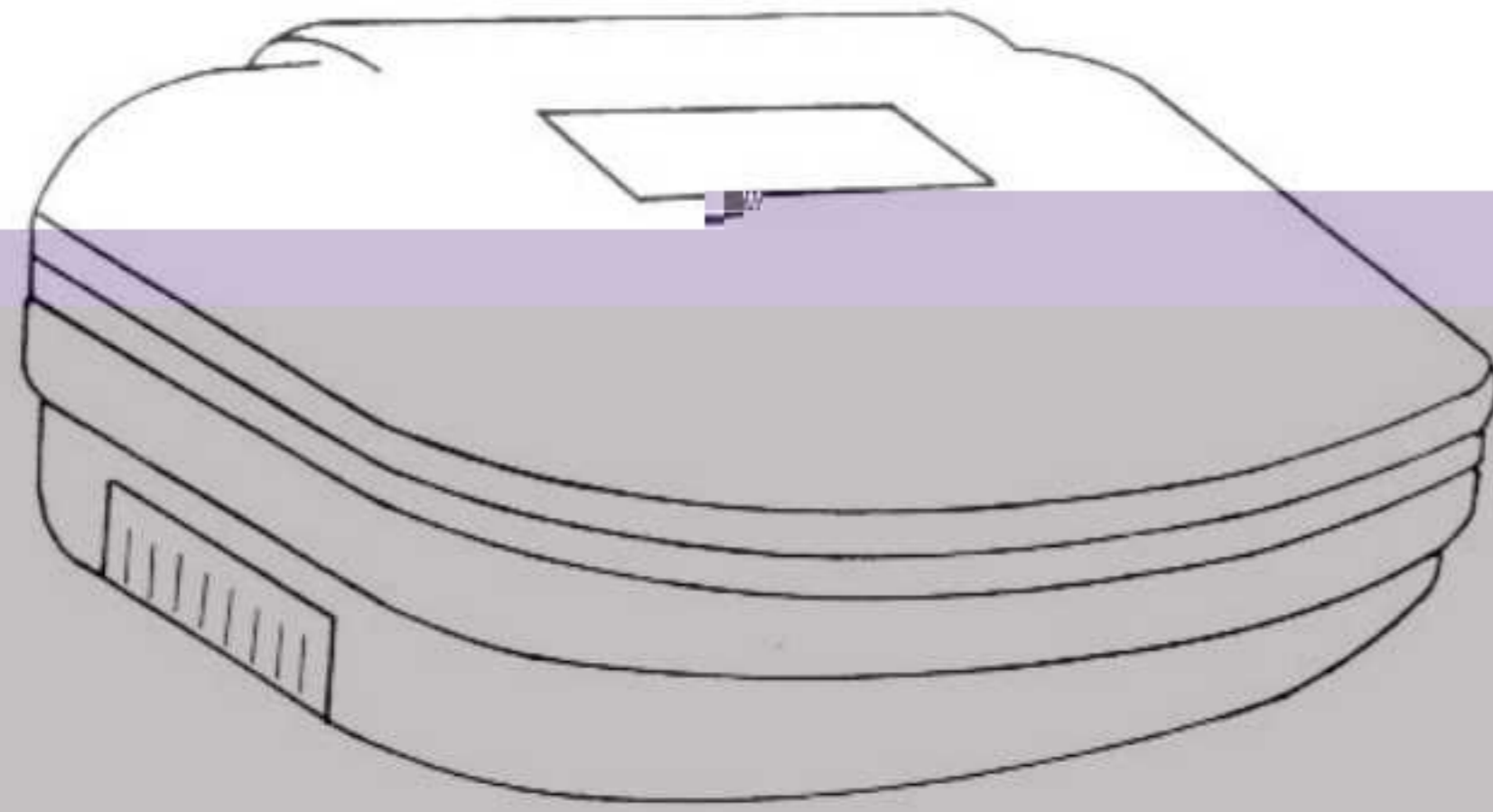
Model J18A₂ -K electric stimulator is a new type one developed by our company based on model J18 for meeting users' requirements. This apparatus optimizes and combines all kinds of frequencies and waveforms, designs 15 multi program mid frequency pulse current prescriptions for massaging various common diseases and frequently occurring diseases, and memorizes them in the processor. In one prescription, all kinds of parameters can be exchanged

several times. Different electric stimulation effects appear. This mid frequency

pulse current stimulation therapy is a comprehensive application of audio

frequency current therapy, pulse modulation mid frequency pulse current therapy and sine modulation mid frequency pulse current therapy, The massage effects are better than mid frequency pulse or low frequency pulse current therapy with a single parameter. You can feel a massage, knock, stir, roll or a quiver when using this apparatus. It is easy to operate and suitable for small hospital and clinic, or in the home.

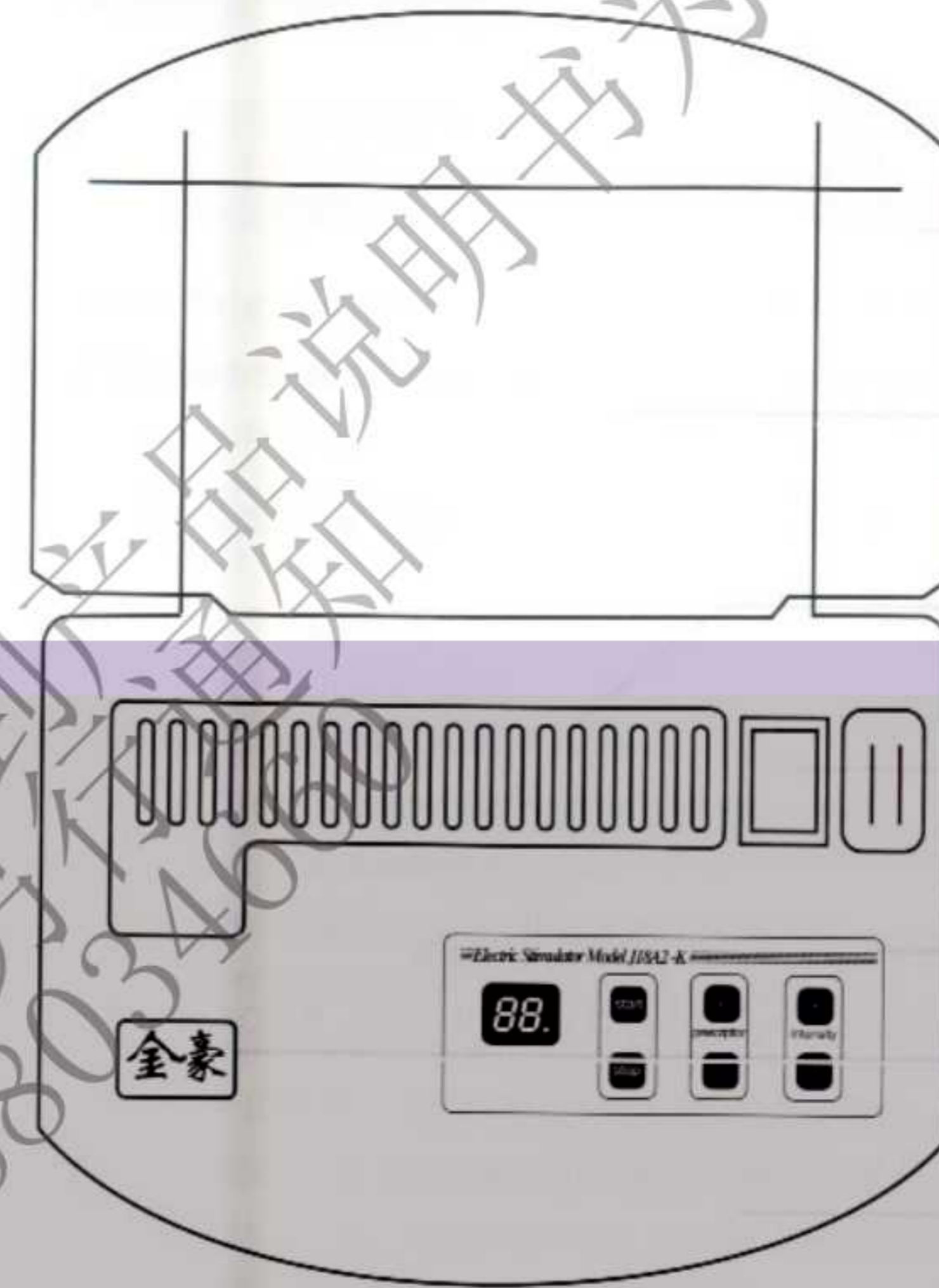upper cover

description

hole for dissipating heat

main supply switch

operating panel

name plate

case

front view of the apparatus

nail for wire winding

socket and fuse

back cover

electrode wire

electrode

back view of the apparatus

II .Operating Panel

1. Nixie Display:

- a. When the indicating lamp does not flash, the prescription number will be indicated
- b. When the indicating lamp flashes, the massage intensity will be indicated.

2. "Start" Key: Press this key, massage starts.

3. "Stop" Key: Press this key, massage stops.

4. Prescription Selection Key: Press "+", prescription number increases progressively; Press "-", it decreases progressively.


III .Operating Step

* Read this operating step in detail and follow it strictly. Do not change sequence with authorization.

1. Open the back cover (see Fig.) and turn on the power supply.
2. Take out the electrode wire from the back cover. One end of wire is connected with two electrode plates.
3. Open the upper cover and push the supply switch. When the apparatus sounds, the nixie display shows No. 1 prescription.


4. Select prescription

Press "+" or "-" key

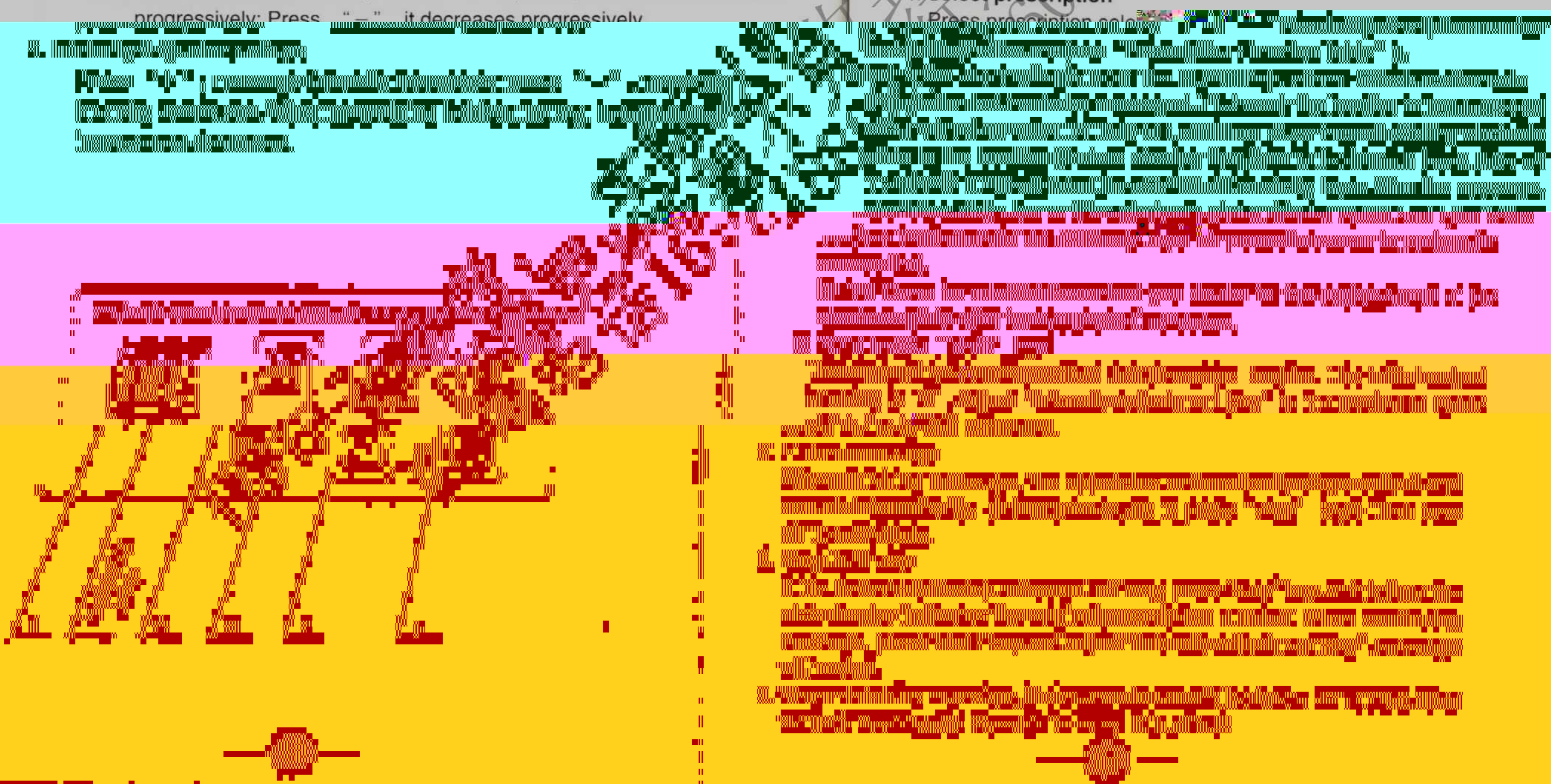※ Attached explanation for the use of remote controller of J18A₂-K model of Electric Stimulator

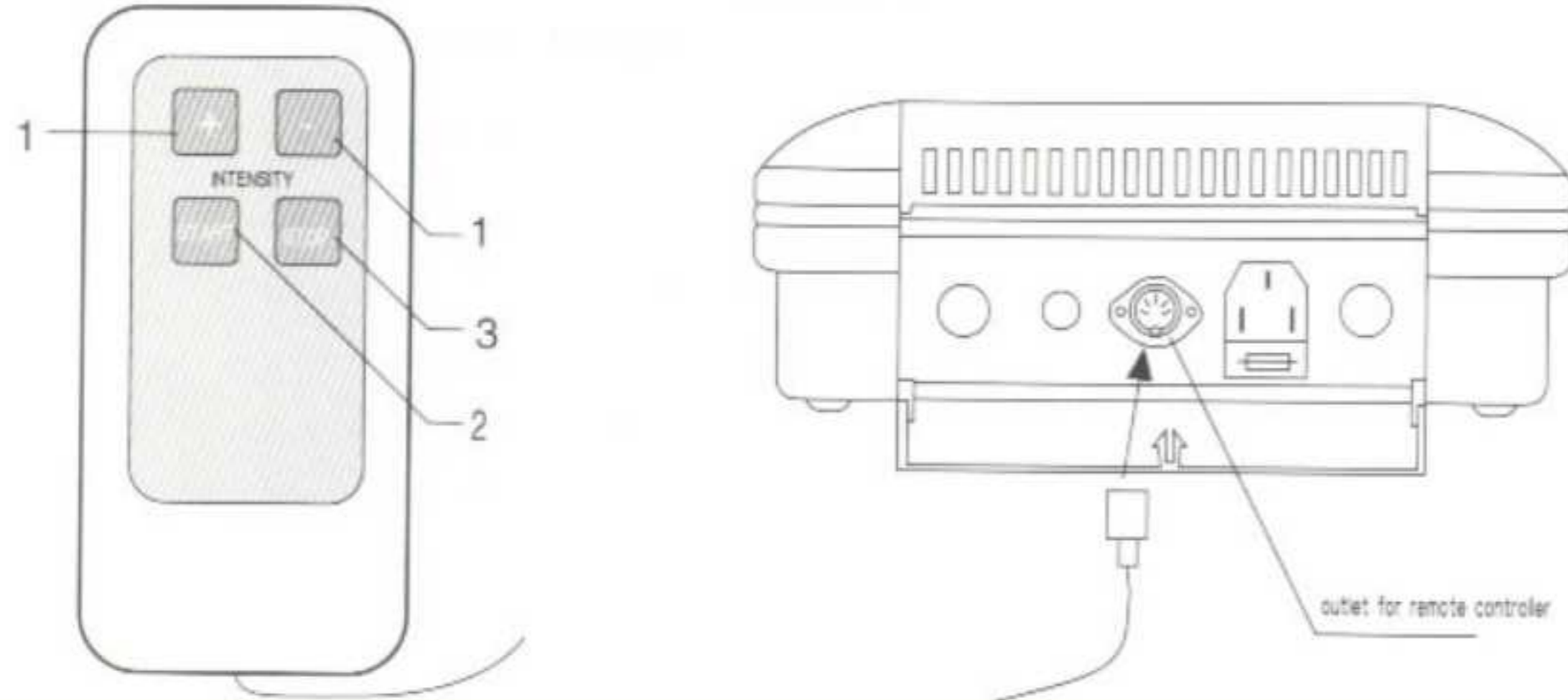	Time	Midfrequency	Low frequency	Modulation waveform
	20min	2kHz	1/2-150Hz	square wave, exponential wave
	20min	4kHz	1/2-150Hz	square wave, exponential wave, sine wave, triangular wave
lylitis.	20min	2kHz	1/2-150Hz	sharp wave, triangular wave, exponential wave
	20min	2kHz	1-150Hz	square wave, exponential wave

IV Prescription Function Instruction

1. Prescription Function Table

Prescription number	Main use	Time	Midfrequency	Low frequency	Modulation waveform
1	Cervical spondylopathy, lumbodynia, sprain, contusion, protrusion of lumbar intervertebral	20min	2kHz	1/2-150Hz	square wave, exponential wave
2	hyperosteoegeny, arthralgia, rheumatoid arthritis.	20min	4kHz	1/2-150Hz	square wave, exponential wave, sine wave, triangular wave
3	scapulohumeral periarthritis, tenosynovitis, external humeral epicondylitis.	20min	2kHz	1/2-150Hz	sharp wave, triangular wave, exponential wave
4	analgesic: used for soft tissue injury	20min	2kHz	1-150Hz	square wave, exponential wave, triangular wave
5	neuritis, sciatica, neuritis of lateral cutaneous nerve of thigh.	20min	3kHz	1/2-150Hz	square wave, exponential wave, sine wave, triangular wave
6	pelvic inflammation, adnexitis, scleroma after injection.	20min	5kHz	1/2-150Hz	sine wave, exponential wave
7	abdominal weight reduction, hip weight reduction.	20min	2kHz	1/3-150Hz	exponential wave, sawtooth wave, square wave
8	functional electric stimulus, hemiplegia.	20min	4kHz	1/5-150Hz	square wave, exponential wave, triangular wave
9	relieving inflammation, subsidence of swelling.	20min	2kHz	1-150Hz	sharp wave, exponential wave, triangular wave

3. Examples of Prescription Application

1 for lumbago

Select No.1 prescription.
Take two 107 × 72mm electrodes
and stick them on back waist
(see Fig.1)

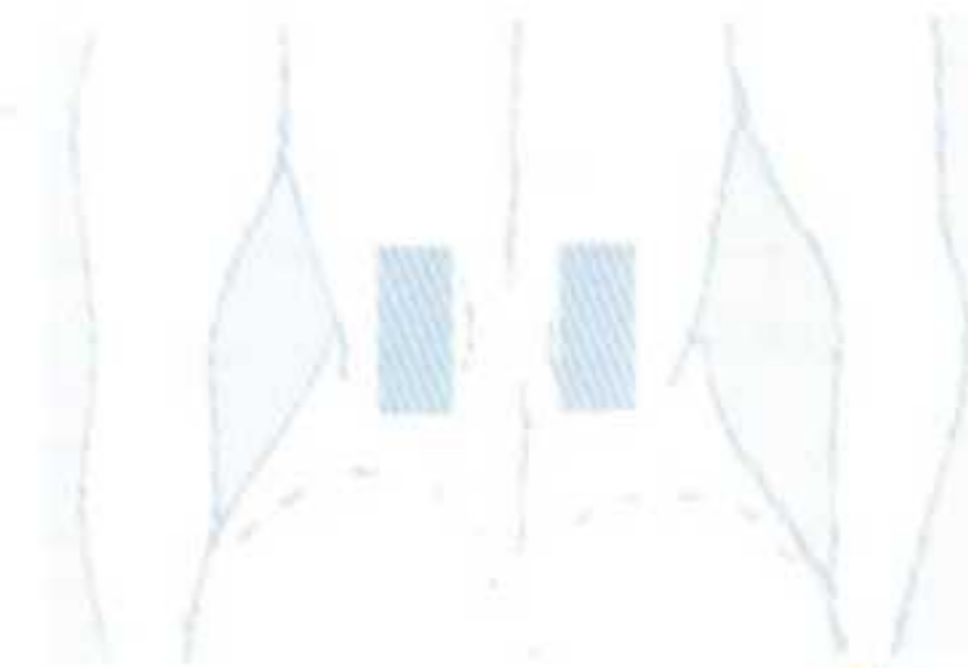
Fig.1

2 for periarthrititis of shoulder

Select No.3 prescription.
Take two ϕ 80mm electrodes
and stick them on front and
back shoulder joints (see Fig.2)

Fig.2

3 for cervical vertebra disease

Select No.1 prescription.
Take two ϕ 40mm electrodes and
stick them on neck.
(see Fig.3)

Fig.3

4 for external epicondylitis of humer

Select No.3 prescription.
Take two ϕ 80mm electrodes and
stick them on internal and
external elbow joints (see Fig.4)

Fig.4

5 for swollen knee joint

Select No.2 prescription
Take two 107 × 72mm electrodes
and stick them on left and right
sides of knee joint. (see Fig.5)

Fig.5

6 for quadriceps m. of thigh atrophy

Select No.12 prescription.
Take two 107 × 72mm electrodes and stick
them on upper and lower ends of
quadriceps m. of thigh of leg (see Fig.6)

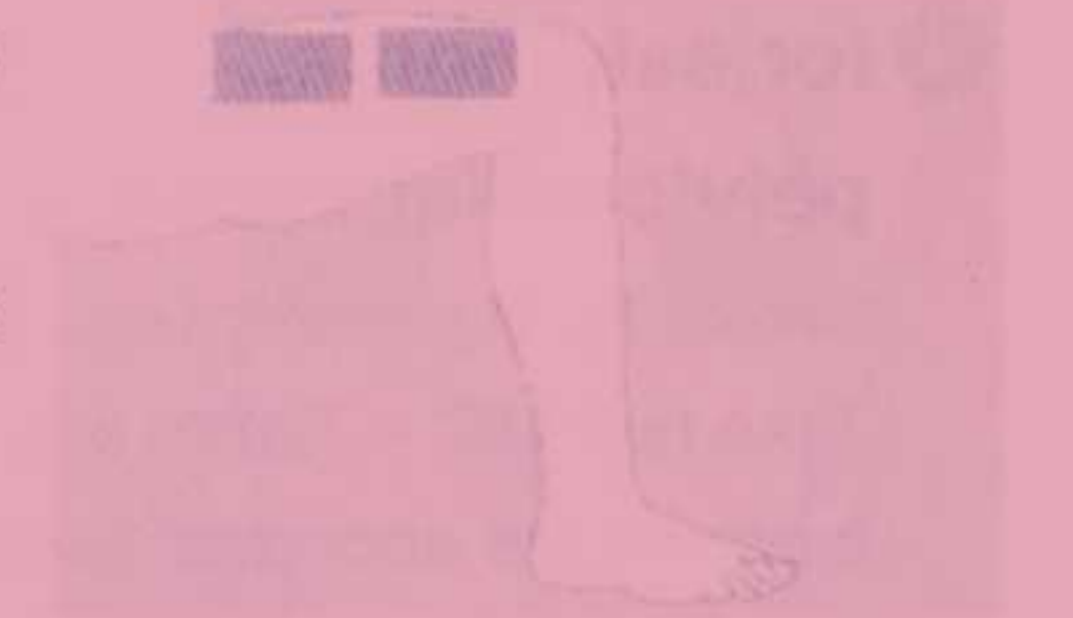
Fig.6

7 for sprain of ankle joint

Select No.1 prescription.
Take two ϕ 40mm electrodes and stick
them on the internal and external
sides of ankle joint (see Fig.7)

Fig.7

8 for hemiplegia (for arm)

Select No.8 prescription.
Take two 107 × 72mm electrodes and stick
them on brachial triceps m. and 1/3 position
above radius (see Fig.8A)

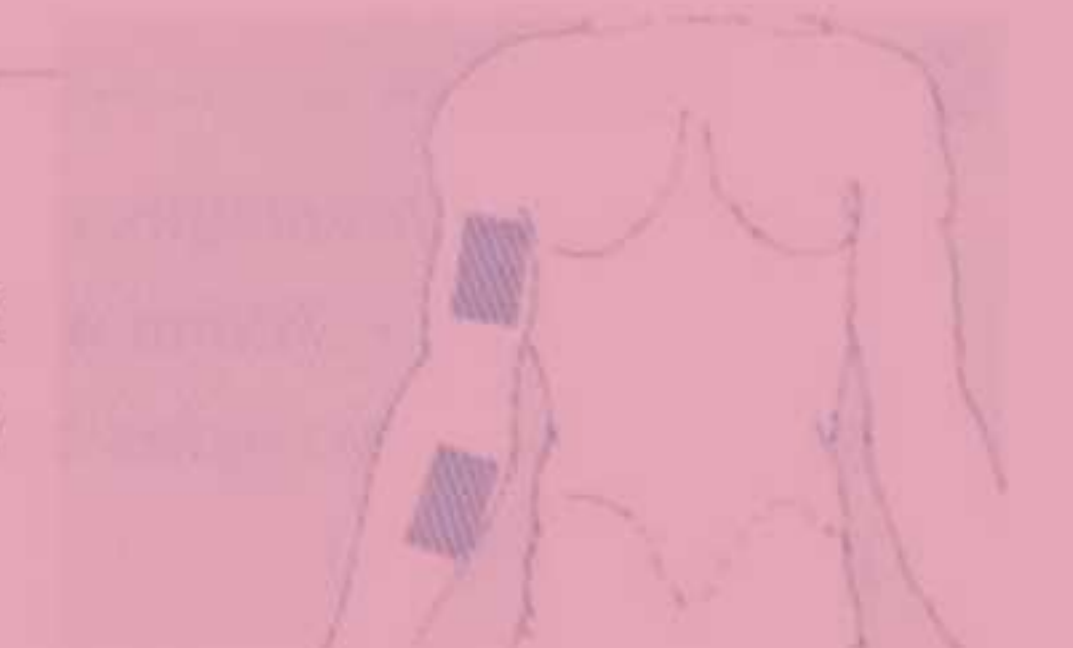
Fig.8A

(for leg)

Select No.8 prescription.

Take two 107 × 72mm electrodes and stick them on quadriceps m. of thigh above patella and front tibia (see Fig.8B)

Fig.8B

9 for salpingitis, pelvic inflammatory disease

Select No.6 prescription.

Take two 107 × 72mm electrodes and stick them on the abdomen (see Fig.9).

Fig.9

10 for gastroptosis

Select No.14 prescription.

Take two 107 × 72mm electrodes and stick them on the belly and the upper of back waist (see Fig.10).

Fig.10

11 for keloid

Select No.15 prescription.

Take two 107 × 72mm electrodes and stick them on two sides of the affected part.(See Fig.11)

Fig.11

12 for laryngitis, vocal cords paralysis

Select No.11 prescription

Take two ϕ 20mm electrodes and stick them on the throat (see Fig.12)

Fig.12

13 for reducing of local obesitas (reducing abdominal obesitas)

Select No.7 prescription.

Take two 107 × 72mm electrodes and stick them on two sides of abdomen (see Fig.13A).

Fig.13A

(reducing of gluteal obesitas)

Select No.7 prescription.

Take two 107 × 72mm electrodes and stick them on the two sides of glutea (see Fig.13B).

Fig.13B

14 for allay, reduction of swelling:

Select No.9 prescription

for pain or germfree inflammation.

Take two proper electrodes and stick them on the pain or swollen position.

15 for sprain, contusion :

Select No.1 prescription.

Stick two electrodes on the sprain
or contusion position.

16 for spur:

Select No.2 prescription.

Stick two proper electrodes
on the spur position.

17 for tendovaginitis:

Select No.3 prescription.

Stick two proper electrodes
on the position of tendovaginitis.

18 for soft tissue injury:

Select No.4 prescription.

Take two proper electrodes
on the injury position.

19 for sciatica:

Select No.5 prescription.

Stick two 107 × 72mm electrodes
on the ischium position
(in the middle of back glutea).

20 for hard mass after injecting:

Select No.6 prescription.

Stick two proper electrodes
on the hard mass position.

21 for eliminating fatigue after sports:

Select No.10 prescription.

Stick two proper electrodes
on the fatigue position
(such as ache position of muscle)

**22 for stimulating and
making muscle developed:**

Select No.12 prescription.

Stick proper electrodes on the local muscle
(such as deltoid muscle of shoulder, bicepsm.
of arm external side muscle of glutea)

23 for peripheral nerve injury:

Select No.13 prescription.

Follow the doctor's advice.

Parallel Penetration massage

Relative Penetration massage

V. Contraindication

Do not use the apparatus under the following conditions: abdomen of pregnant woman, cancer, heart, tuberculosis, acute pyogenic inflammation, hemorrhage position, thrombophlebitis, tetanus, larger metal foreign matter in the massage position and one wearing cardiac pacemaker.

VI. Caution


1. When placing the electrode plate on the human body, do not turn the supply on or off.

2. When the apparatus is working, it must not be near the refrigerator, washing machine, microwave oven and other high frequency devices. To prevent disturbance and over intensive electric stimulation, do not share one socket with these devices.

3. The self-sticking electrode plate can neither be sterilized by boiling or by alcohol. Wash it with clear water after each massage, shake off the droplets and stick the plastic membrane on it. It is advisable to keep it only for personal use. As the viscosity and conductivity of the electrode plate shall become depressed after repeated use, new electrode pieces shall be replaced, and it is advisable to do so after it being used for 20 times.

4. Two ends of output are not allowed to come into contact during therapy because short circuiting and damage of the apparatus will occur.

5. Dry skin and secretion can cause poor contact and affect massage effects. So

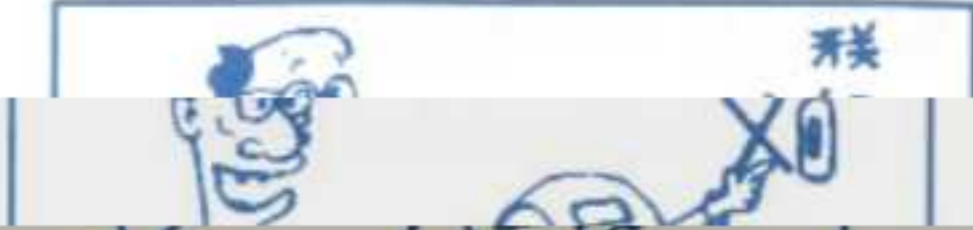skin should be cleaned and kept wet before massage.

6. In the massage course strong or weak or even no current is normal and is set by prescription.

7. If after turning on the supply, the prescription number and indicating lamps are not lighted, first check the supply socket and fuse.

8. When AC voltage is bellow 190(or 90V) or over 240(or 130V), the apparatus will not normally operate.

9. The electrode plate must not come into contact with oil products.

10. When there is acute sprain, do not use the apparatus until 24 hours after the sprain occured.

3. output waveform: Sine wave, square wave, triangular wave, exponential wave, sharp wave, saw tooth wave, combined wave.
4. percentage modulation: 100%, 75%, 50%, 25%, 0%.
5. output peak current: $50\text{mA} < \text{current} < 100\text{mA}$ (load 500Ω)
6. modulation mode: continuous, intermittent, variable
7. intensity adjustment: continuous
8. working voltage: AC 220V (or 110) $\pm 10\%$, 50/60Hz
9. power: $\leq 15\text{W}$
10. output channel: single

Attachments Table with the Apparatus

Number	Name	Quantity
2	main apparatus	1
3	supply wire	1
4	electrode plate ① 107 × 72 ② φ 80 ③ φ 40 ④ φ 20	4pairs
5	bandage	1
6	manual	1
7	certificate	1
8	guarantee card	1

VI. Main Technical Parameters

1. mid frequency pulse frequency: 1-10KHz
2. low frequency modulation frequency: 1/8-150Hz